



eTwo®

Complete Face &
Body Renewal



Complete Skin
Renewal



elōs
results



elōs
safety

"In my trials, I have found that the Sublative treatment from Candela offers a unique combination of ablation and skin resurfacing leading to wrinkle reduction with less downtime than traditional skin rejuvenating treatments.

My patients are pleased with the fast, visible results of smoother and brighter skin, without having to sacrifice time from their busy schedules."

Amy Taub, M.D.
Dermatologist, Illinois

Treat acne scars & striae

Treat facial wrinkles

eTwo combines the unique Sublime and Sublative applications for comprehensive facial rejuvenation in as little as a single visit.

Sublime uses the eIōs combination of infrared light and bi-polar radiofrequency energies to remodel the deep dermal layer of the skin by stimulating new collagen growth.

- Improvement of wrinkles
- Smoother skin appearance
- Based on Candela's eIōs technology - a combination of energies that stimulates new collagen growth¹

Sublative fractionated bi-polar radiofrequency technology places the heat energy successfully into the upper dermis where it can produce a significant increase in collagen with minimal epidermal disruption.

- Improvement of acne scars and striae
- Treat facial wrinkles
- Use on ALL skin types

Sublime and Sublative Complementary Therapy

Wrinkles



Post 2 treatments
Photos courtesy of Ruthie Amir, M.D.

Wrinkles and Contouring



Post 1 treatment
Photos courtesy of Brian Biesman, M.D.

Sublime

Wrinkles



Post 2 treatments
Photos courtesy of Ruthie Amir, M.D.


Sublative

Striae



Post 3 treatments
Photos courtesy of Tsai, M.D.

Wrinkles

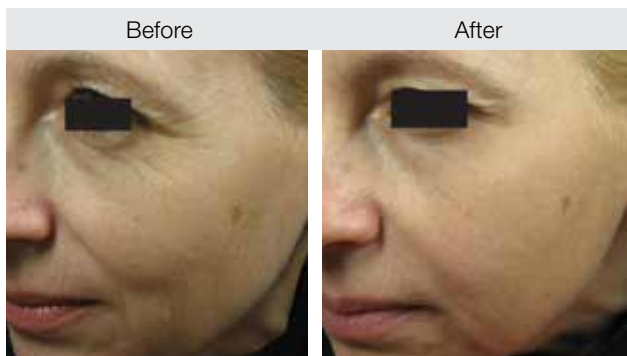

Post 2 treatments
Photos courtesy of Ruthie Amir, M.D.

Acne Scars



Post 3 treatments
Photos courtesy of Sasha Jazayeri, M.D.

"In my experience, when treating patients with eTwo the results I see in skin quality and luminosity are consistently impressive."

Brian Biesman, M.D.
F.A.C.S., Tennessee

Photos have been unretouched. Individual results may vary.



eTwo

Visible Improvement in as Little as a Single Visit

Now eTwo has enhanced energy and a Sublative and Sublime complementary treatment protocol for complete rejuvenation in fewer sessions.

User interface intuitive, easy-to-learn

Easily portable

2 applicators connected at once for easy treatment

Touchscreen 15 inch and foldable



Sublime Applicator



Sublative RF Applicator

The Technology of Sublime and Sublative

Sublime & Sublative: The Versatile Advantage

Sublime uses the eIōs combination of IR and Bi-Polar RF to provide controlled deep dermal heating. This predictable heating of the deep dermis both stimulates and shortens aged collagen fibers, resulting in the improvement of wrinkles and facial contour.

The unique Sublative ergonomic applicator delivers bi-polar radio frequency energy via the matrix of electrodes of the disposable Sublative iD tip. This fractional delivery of RF energy causes micro-epidermal ablative injuries with wider spread dermal coagulation. This initiates a strong wound healing response in the dermal matrix while also minimizing the damage to the epidermis. Less affected tissue surrounding the matrix points may be stimulated by heat which helps accelerate the healing process, leading to a smoother, rejuvenated appearance of the skin.



Su-blime
[suh-blahym], noun.


1. Non-invasive wrinkle treatment
2. Uses Candela's eIōs technology
3. Fast treatment with visible results



Sub-la-tive
[sah-bla-tiv], noun.

1. A new category of fractional skin treatments
2. High dermal impact, low epidermal disruption for all skin types
3. Minimal downtime, proven success, satisfied patients


The Science of Sublime

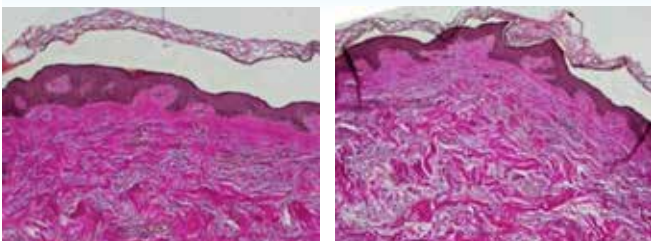

Pre treatment

Deep dermal remodeling 3 months post treatment

The Science of Sublative

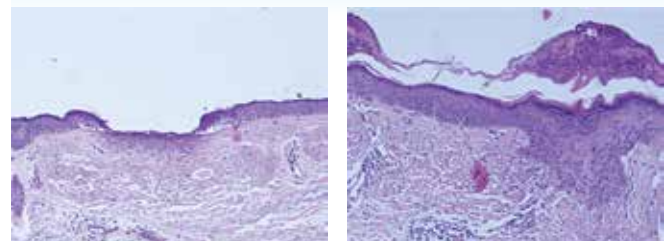

Histology showing single electrode immediate Sublative impact

Histology showing rapid healing and strong dermal remodeling only 4 days post a single treatment



Visible Results, Low Downtime

The promise of noticeable results with minimal downtime and lower procedure cost is spurring interest in fractional methods of skin resurfacing. Patients are increasingly requesting these treatments to improve their appearance.

The eTwo system not only enables your practice to address this growing demand, but also adds non-invasive wrinkle treatment to your revenue stream.

Meet the Growing Demand

The eTwo system was designed with your practice in mind. Excellent ROI potential - lower acquisition and disposable costs, plus lowest overall cost of ownership compared to competitive products.

eTwo System – A Smart Choice

1. Intelligent Feedback System provides real time information on skin condition, impedance and the energy delivery of each pulse for consistency and reproducibility.
2. Minimal downtime so patients can quickly resume their normal activities.
3. Suitable for all skin types.¹
4. Quick treatments: typically up to 1 hour depending on area treated with a Sublative and Sublime complementary protocol.
5. Short treatment protocol: Visible improvement after the first session.
6. Lightweight, compact device that easily transports between locations.



Sublative RF™

Multiple Treatment Options, Exceptional Results - Dramatically improve skin quality, tone, wrinkles, acne scars and striae.¹

Sublative iD Tips

1. Intelligently designed tips to maximize treatment options and results
2. Re-connectable per treatment
3. Multiple tip configurations



Sublative iD Standard 64 pin tip
100 and 400 pulses per tip



Sublative iD Focal 44 pin tip
200 pulses per tip


64 pin

1. Fast full face Sublative resurfacing for improved skin tone and quality
2. Multiple tip configurations for various treatment protocols

44 pin

1. Optimized electrode geometry for smaller areas and focused treatment on peri-orbital and peri-oral wrinkles
2. Focused treatment of acne scars and forehead wrinkles

Sublative RF Treatment Screen

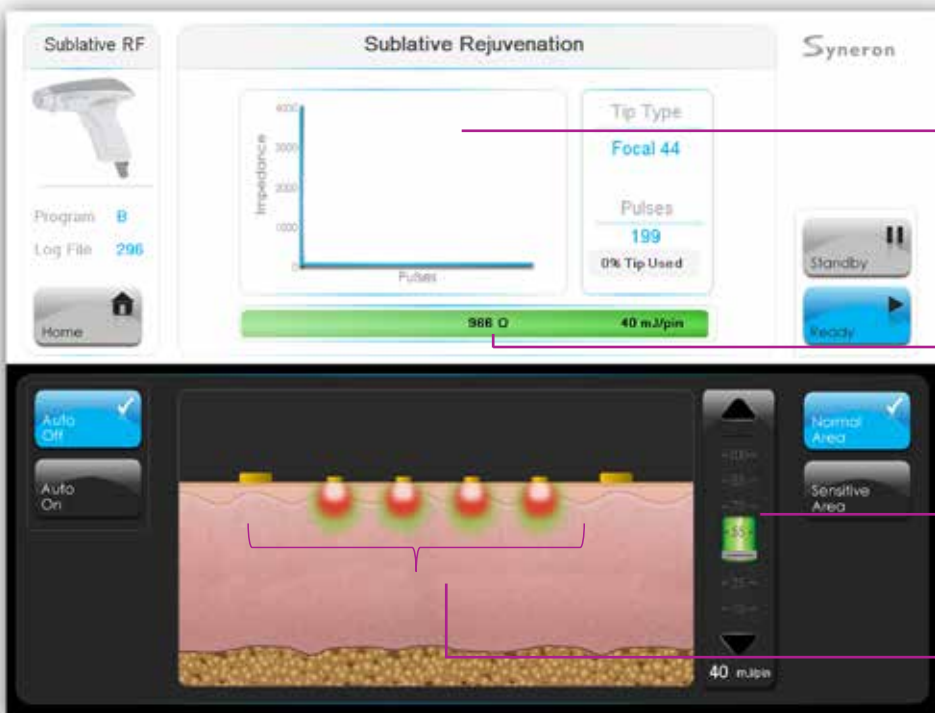

Intelligent Feedback Graph
Shows a moving average of pulse to pulse impedance

Skin Feedback Meter
Shows skin condition and energy delivered

Full Energy Selection Bar
With highlighted program range

Scaled Impact Pattern
Shows depth and width of thermal damage



System Specifications

Sublative RF Energy	Up to 100 mJ/pin
Sublative RF iD Spot Size	Standard 64 - 12 x 12 mm Focal 44 - 11 x 3 mm
Disposable Sublative RF iD Tips	Standard 64 up to 400 pulses Focal 44 up to 200 pulses
Sublime Energy	Up to 200 J/cm ³
Sublime Spot Size	12 x 8 mm
User Interface	15" touch screen
System Dimensions Table Top	53 x 30 x 40 cm 21.8" x 11.8" x 15.7"
System Weight	8.5 kg / 18.7 lbs
Electrical Requirements	100-240 VAC; 3A; 50-60 Hz

eTwo



For more information about how the Two may help achieve your practice goals, contact your local Candela sales professional, or visit candelamedical.com.



Disclaimer: All contents of this material are for informational purposes only and provided by Candela without warranties of any kind. Healthcare professionals are solely responsible for making their own independent evaluation as to the suitability of any product for any particular purpose and in accordance with country specific regulations. The availability of products and the indications mentioned in this material is subject to the regulatory requirements and product registration status in each country. Refer to the User Manual for country specific indications. Products and technical specifications may change without notice. Please contact Candela for more details.

CE 0344

© 2019 Candela Corporation. This material contains registered and unregistered trademarks, trade-names, service marks and brand names of Candela Corporation and its affiliates. All other trademarks are the property of their respective owners. All rights reserved. PB75230EN